



Special Edition


AccuRay 525R

Innovative & Economic System for Multi-Purpose X-ray Imaging Solution

DK⁺ Medical Systems

AccuRay-525R


World's highest specification High Grade inverter technology. DK has produced an X-ray high voltage Generator with a high frequency Inverter. This High Frequency Inverter produces low ripple tube voltage and high X-ray quantum efficiency. The frequency is variable according to load conditions. A stable X-ray output is achieved even with increasing tube currents.

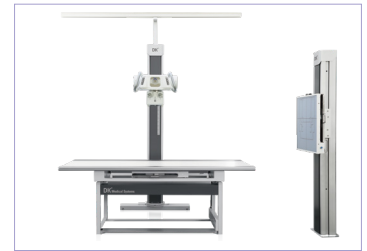

|| AccuRay 525R

Specifications

System	AccuRay 525R
Radiography kV	40~125kV
Radiography mA	10~500mA
Generator Time	0.004~10 sec
Self-diagnostic function	LED error code display
Anatomical Programming	216 programs
Disply System	Digital LED display
Automatic Exposure Control (AEC)	Optional

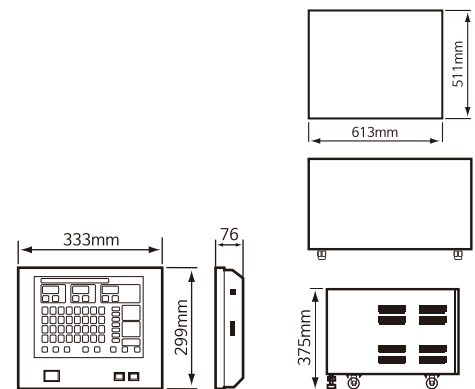


Console size

Generator size

Structure



|| SF-80B

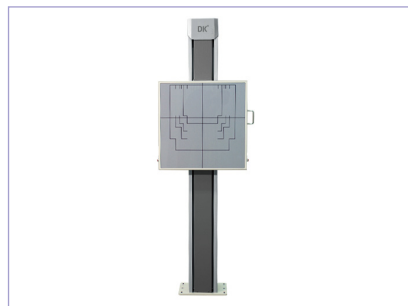

|| VB-20



|| TF-32N

